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**GB – Tip:** The dispensing chemist you trust can help you with the setting up of the mediTimer<sup>®</sup>, e.g. how to handle half tablets, medication that needs to be taken a number of times per day but only exists once or medication that is to be stored in the refrigerator. You will find additional information at [www.mediTimer.de](http://www.mediTimer.de).

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**GB – Thank you for choosing the mediTimer®. You can now organise your health safely!**

The mediTimer® is a simple and intelligent organisation system that supports you with the taking of your medication and safely leads you through the day.

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**GB – The System at a Glance**

**mediTimer®**

The basic module with  
3 medication boxes,  
2 sheets of labels with the  
administration times

**Available separately:**

the expansion module with  
1 medication box

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**GB – The benefits of your health manager**

- ✓ 1 x sorted – weeks of peace and quiet!
- ✓ Customised correct times to take the tablets
- ✓ Only safe original packaging
- ✓ For all forms of administration
- ✓ Expandable for additional administration times

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**GB – Available separately**

The mediTimer® can be separately expanded for more than three administration times.

In this case, place the expansion module with the connecting tab of the inserted box facing upward. Now press the mediTimer® onto the expansion module until it clicks.

**Please note:**

the expansion module can no longer be removed!

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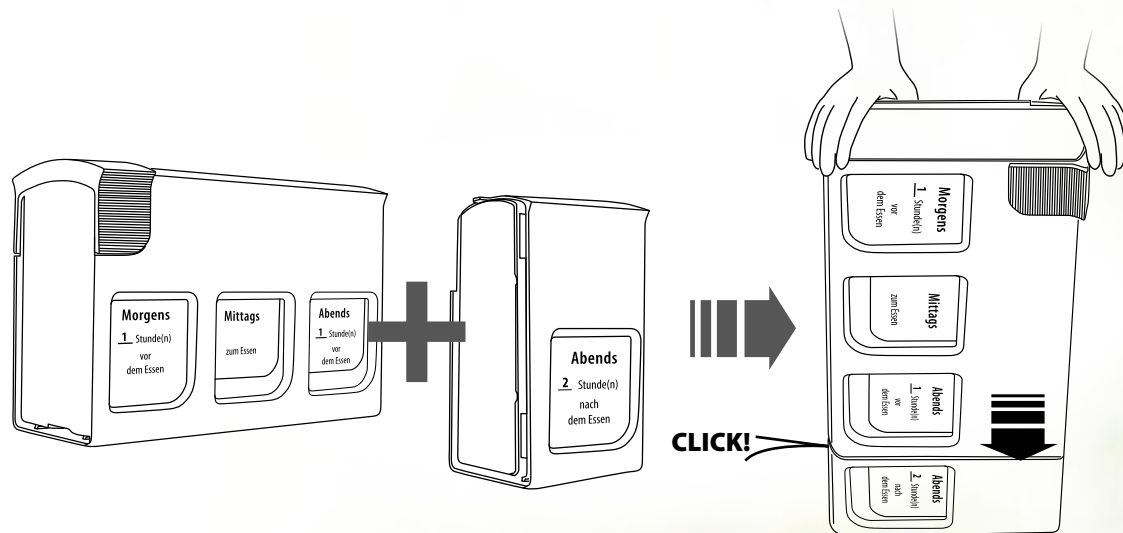
#### GB – Preparation

Select the labels in keeping with the customised correct administration times (e.g. on the basis of your medication regimen), inscribe them and adhere to the boxes.

**Tip:** Should you have any questions about setting the mediTimer<sup>®</sup> up, you can obtain information from our website: [www.meditimer.de](http://www.meditimer.de) or contact us.

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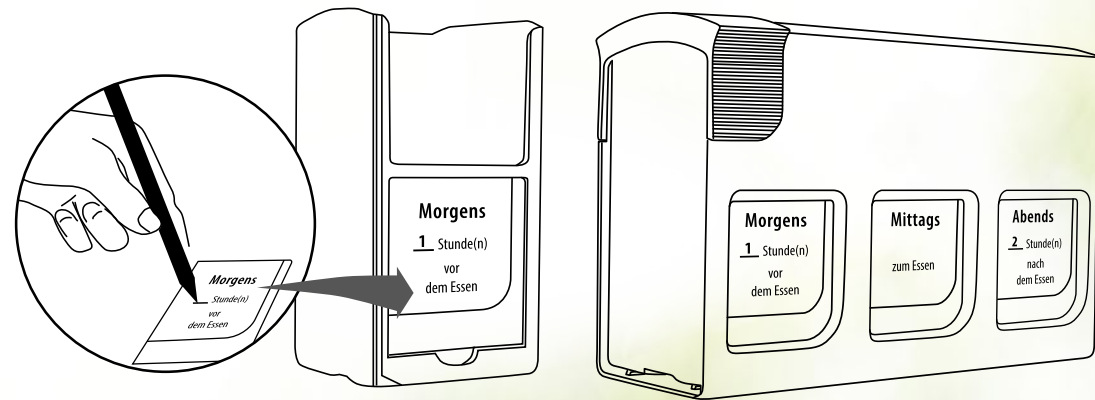
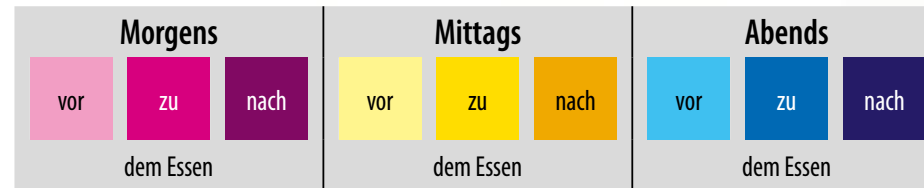
#### GB – Filling the boxes


- 1 Pull the barrier up.
- 2 Push all of the boxes out from right to left.
- 3 Fill the boxes with your medication in keeping with the appropriate times. Everything that is to be taken/applied at a certain time is to be placed in the correct box once. For example, one tablet strip (original blister) for tablets or the drops bottle with the label stating how many drops are to be taken for drops. Should there only be one blister for two administration times, this is to be cut through and distributed between the two appropriate boxes.  
**Tip:** Should you have any questions about filling the mediTimer® up, you can obtain information from our website: [www.meditimer.de](http://www.meditimer.de) or contact us.

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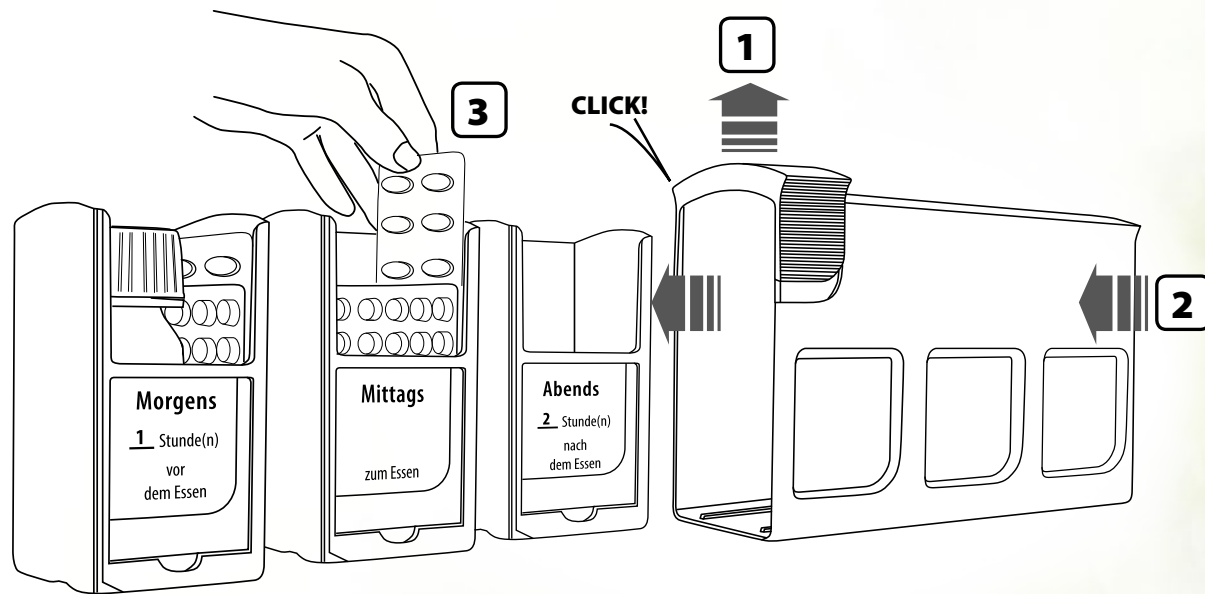**GB – Filling the mediTimer® with the filled boxes**


- 1 Push the barred down.
- 2 Insert the filled medication boxes in the housing from right to left in the correct day order.

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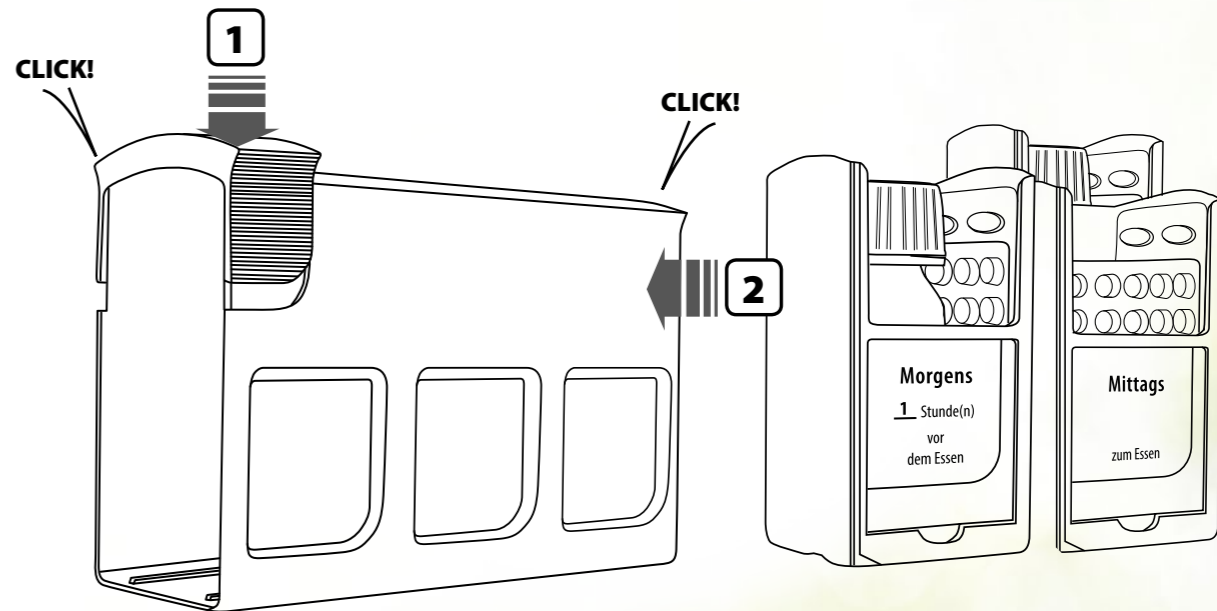#### GB – Use

- 1 Pull the barrier up.
- 2 Push the first box to the front at the required times.
- 3 Take/ apply the medication.
- 4 Sperre wieder runterdrücken.
- 5 Push the barrier back down.
- 6 The box with the next administration time is now automatically at the front.

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