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Information
Information
Information
Informazione

Morgens

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Mittags

zum Essen

Abends

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GB – Many thanks for your interest in the mediTimer®. You can now organise your health with the mediTimer®.

The mediTimer® is a simple and intelligent organisation system that supports you with the taking of your medication and safely leads you through the day. On the next pages, you will be able to comprehensively familiarise yourself with the mediTimer as a completely new generation of medication dosing systems.

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GB – The System at a Glance

mediTimer®

The basic module with
3 medication boxes,
2 sheets of labels with the
administration times

Available separately:

the expansion module with
1 medication box

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GB – The benefits of your health manager

- ✓ 1 x sorted – weeks of peace and quiet!
- ✓ Customised correct times to take the tablets
- ✓ Only safe original packaging
- ✓ For all forms of administration
- ✓ Expandable for additional administration times

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GB – Available separately

The mediTimer® can be separately expanded for more than three administration times.

In this case, place the expansion module with the connecting tab of the inserted box facing upward. Now press the mediTimer® onto the expansion module until it clicks.

Please note:

the expansion module can no longer be removed!

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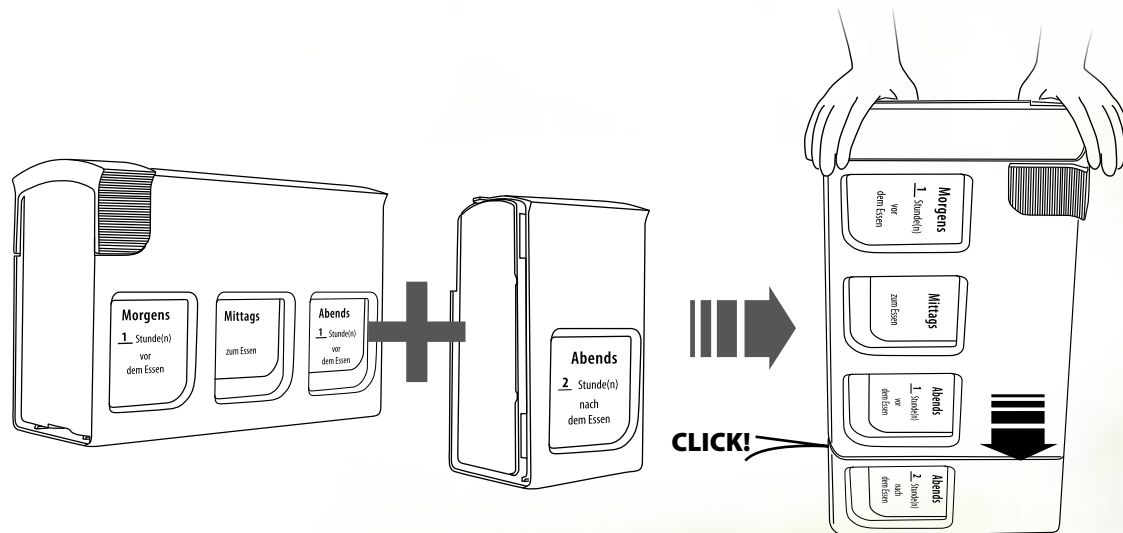
GB – Preparation

Select the labels in keeping with the customised correct administration times (e.g. on the basis of your medication regimen), inscribe them and adhere to the boxes.

Tip: Should you have any questions about setting the mediTimer[®] up, you can obtain information from our website: www.meditimer.de or contact us.

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
EN – Innovative colour guidance system for the correct administration times (further information is available at www.mediTimer.de)

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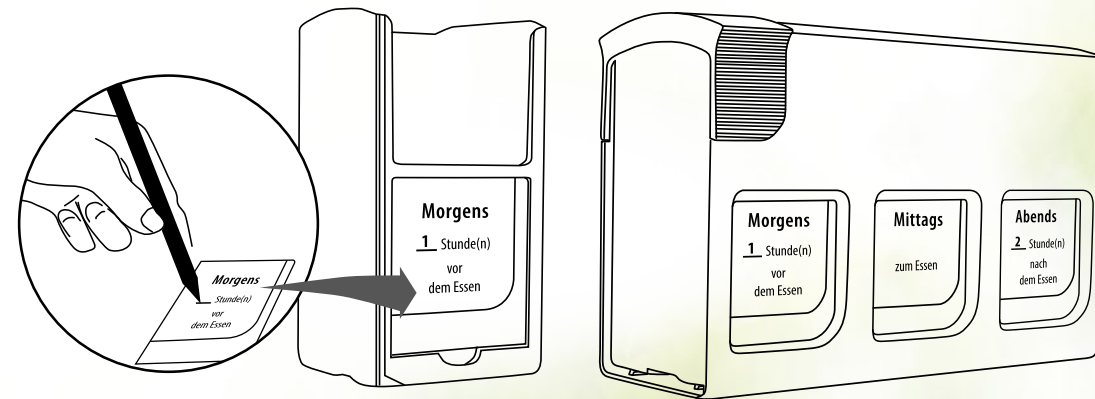GB – Filling the boxes


- 1 Pull the barrier up.
- 2 Push all of the boxes out from right to left.
- 3 Fill the boxes with your medication in keeping with the appropriate times. Everything that is to be taken/applied at a certain time is to be placed in the correct box once. For example, one tablet strip (original blister) for tablets or the drops bottle with the label stating how many drops are to be taken for drops. Should there only be one blister for two administration times, this is to be cut through and distributed between the two appropriate boxes. **Tip:** Should you have any questions about filling the mediTimer® up, you can obtain information from our website: www.meditimer.de or contact us.

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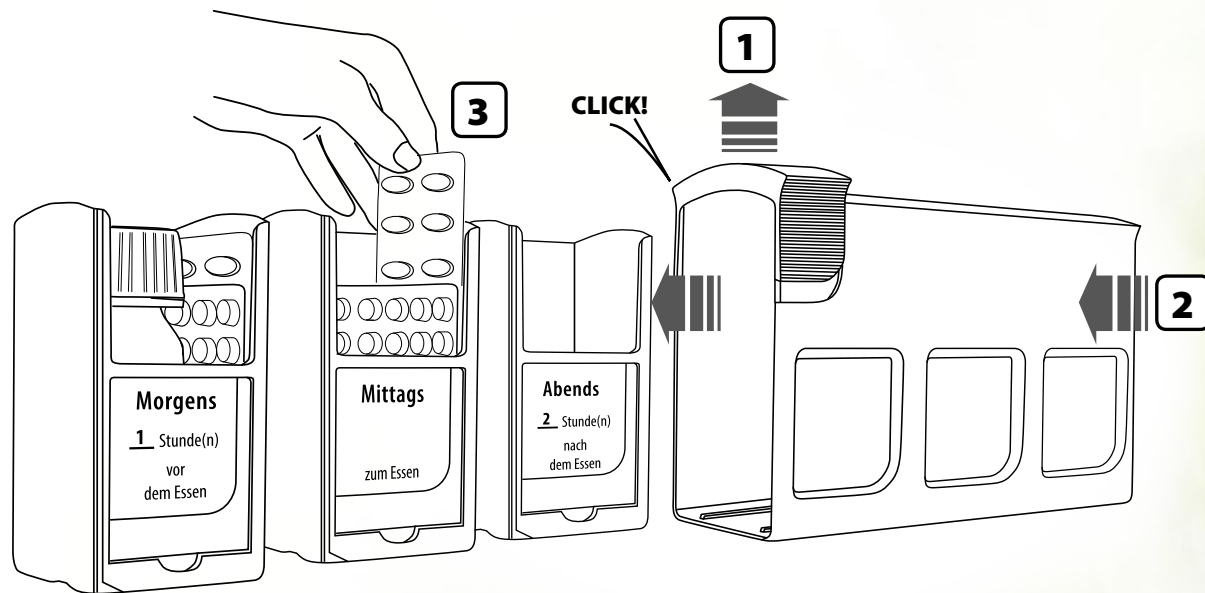GB – Filling the mediTimer® with the filled boxes


- 1 Push the barred down.
- 2 Insert the filled medication boxes in the housing from right to left in the correct day order.

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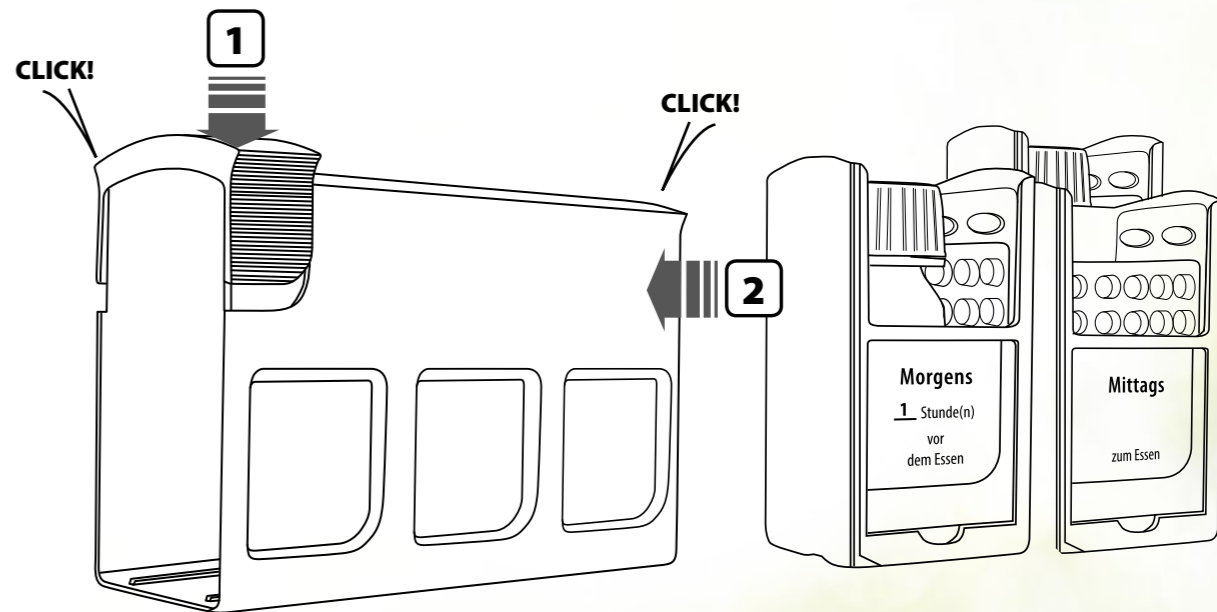GB – Use

- 1 Pull the barrier up.
- 2 Push the first box to the front at the required times.
- 3 Take/ apply the medication.
- 4 Sperre wieder runterdrücken.
- 5 Push the barrier back down.
- 6 The box with the next administration time is now automatically at the front.

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